



Exclusively dedicated  
to the medical field  
since 2002



**STORAGE CABINETS**  
FLEXIBLE ENDOSCOPES & ULTRASOUND PROBES



# Decrease reprocessing frequency and maximize your medical device availability thanks to our storage cabinets

**Allows up to 7 days of storage before reprocessing**

- 100 % stainless steel
- LED Lighting system
- Movable endoscope hooks/supports
- Removable drip tray
- Easily combine large and small diameter endoscopes in 1 cabinet
- Pre-filter and Hepa filter
- Removable basket for caps and valves
- Compliance with Canadian and American health care standards
- Air-tight glass doors with key-lock
- Full traceability system (Optional)
- Security system with lock and integrated alarm (Optional)
- Channel Purge Storage for Faster Drying (Optional)



Made entirely of Stainless Steel.  
Design and Produced in Canada.



The ergonomics and designs of our storage cabinets are the result of collaborative work accomplished with over one hundred users in a real work environment.



Fully customizable hook system. Mixed storage of medical devices of various diameters from all manufacturers.

# Intelligent storage that offers control and traceability of your medical devices



ChronoMEDIC 

Our pre-integrated Traceability of Endoscopes™ and ChronoMEDIC Cabinet Management™ modules are efficient tools which simplify the day-to-day management tasks of your department heads.

In addition to ensuring optimal rotation of your medical devices, they will allow you to decrease your repair costs from overuse of your equipment. Built-in controls and security features help prevent the loss and theft of your flexible endoscopes and ultrasound probes.

- Pre-cleaning, cleaning and disinfection
- Quality control
- Periodic maintenance
- Validation of reprocessing timeline
- Repair shipment and return
- Suggested use order
- General history consultation



## ChronoMEDIC™ management module for storage cabinets


- Complete and secured identification of users
- Identification for control of entry and exit of each endoscope
- Alert and simple identification of endoscopes needing to be reprocessed
- Quick identification of the location of the next endoscope to be used
- A precise and detailed list of endoscopes by storage location
- Complete traceability of endoscope movements to other cabinets
- Consultation of endoscopes by type and by expiry date
- Consultation of complete history of entries and exits by day and by user
- Printing of multiple reports as needed

## ChronoMEDIC™ management module for endoscopes

- Intuitive Control and guidance
- Quality control for each reprocessing step
- Periodic maintenance of endoscope reprocessors according to an established schedule
- Suggestion of which endoscope should be used next (by type, date and hour)
- Pre-cleaning, cleaning and high-level disinfection of endoscopes with multimedia instructions
- Sending and returning of endoscopes for repair with quality control at return, and cost management
- Research of endoscopes used by patients
- Recall management by endoscope or by patient
- Validation of reprocessing time and quality when using for a patient
- Detailed history of endoscope reprocessor usage
- Consultation of the general history (usage, cleaning, washing and disinfection machine, repairs, etc.)

# Models of cabinets regularly used in the department of:

- Gastroenterology
- Respiratory Care
- Urology
- Anesthesia
- ENT Clinic
- Emergency & Operating Rooms
- Ultrasound
- Medical Imaging

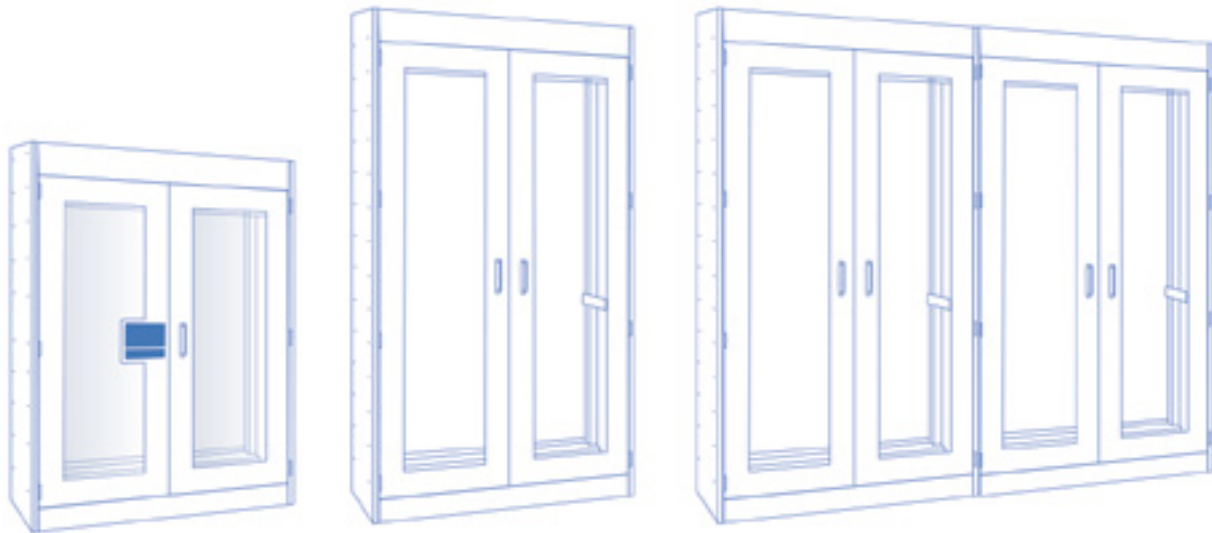

## STAINLESS STEEL STORAGE CABINETS FOR SMALL DIAMETER ENDOSCOPES

bronchoscopes, nasolaryngoscopes, cystoscopes...

- Model 5 endoscopes - (H-62 3/4" x L-29 3/4" x D-18")
- Model 10 endoscopes - (H-62 3/4" x L-49 1/2" x D-18")
- Model 12 endoscopes - (H-62 3/4" x L-35 1/4" x D-18")
- Model 18 endoscopes - (H-62 3/4" x L-49 1/2" x D-18")

## STAINLESS STEEL STORAGE CABINETS FOR LARGE DIAMETER ENDOSCOPES

colonoscopes, gastroscopes, duodénoscopes...

- Model 5 endoscopes - (H-92 1/4" x L-29 3/4" x D-18")
- Model 10 endoscopes - (H-92 1/4" x L-49 1/2" x D-18")
- Model 12 endoscopes - (H-92 1/4" x L-35 1/4" x D-18")
- Model 18 endoscopes - (H-92 1/4" x L-49 1/2" x D-18")



We also manufacture custom cabinets according to your specific dimensions and configurations



625 – 1 Av. Holiday, Tour Est, Pointe-Claire, Quebec, Canada H9R 5N3  
Tel.: 514-697-8087 | 1-877-229-2224 Fax: 514-697-0993  
info@focusendoscopy.com